

LUMINETTE Q

The perfect addition for highly effective laser tattoo removal treatments

British Manufactured Q-Switched Nd:YAG Laser at 532nm & 1064nm



Lynton

MADE IN
BRITAIN

Unrivalled Passive Q-Switched Laser for Professional Tattoo Removal



From the very start, Lynton has been a pioneer in the science and manufacture of laser tattoo removal devices. As a result, we have encapsulated our vast experience and expertise into one, cost-effective system - the LUMINETTE Q.

Purposely designed to fit straight into any tattoo studio, clinic or salon, this portable, desktop-based passive Q-switched Nd:YAG laser will help you achieve unrivalled clearance and reduction of professional, amateur or traumatic tattoos.

High-quality British manufacturing, together with our award-winning support network (training, technical, clinical & business), ensures the LUMINETTE Q will become a highly-profitable asset for your business.

Over 25 years experience in the development & manufacture of world-leading tattoo removal technologies

Intelligent laser handpieces with intuitive, easy-to-use software to provide full control over treatment parameters

Cost-effective passive tattoo removal laser with an unrivalled clearance of tattoos, together with **award-winning post-purchase support & service**



LUMINETTE Q interchangeable laser spot sizes




The 1064nm handpiece for treatment of black & dark blue ink pigments

WINNER 2017
Aesthetics Awards
 EQUIPMENT SUPPLIER OF THE YEAR



The Science?

The LUMINETTE Q generates high-powered, extremely short (nanosecond) pulses of laser light which are then absorbed by the tattoo ink particles, breaking them down into smaller pieces. These shattered pieces are then absorbed by the body's immune system, which disperses the ink. A gradual fading of the tattoo occurs over a series of treatments.

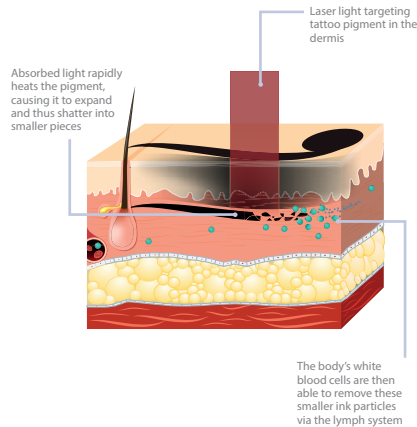

Customer Review



Since working with the Lynton LUMINETTE Q system, we haven't looked back! We're getting much better results compared to our previous laser and much more business because of Lynton's reputation and the results their systems provide. The laser is easy to use and the training was extremely thorough. We find any enquiries are dealt with by Lynton very efficiently. We would recommend Lynton to anyone. We certainly wouldn't ever be using another machine.



Results

Tattoo Removal



Before

After

Tattoo Removal

After 5x treatments



Before

After

Tattoo Removal

After 4x treatments



Before

After

Tattoo Removal

After 3x treatments



Before

After

Tattoo Removal

After 4x treatments



Before

After



Call Lynton today on 01477 536 977 or
visit lynton.co.uk

System Specifications LUMINETTE Q

Laser Type	Q-Switched Nd:YAG	
Wavelength	1064nm	532nm
Pulse Width (ns)	6	6
Repetition Rate (Hz)	Up to 10	Up to 10
Spot Size (mm)	3 & 5	3 & 5
Fluence Range (J/cm ²)	1.2 – 16.9	1.2 – 16.9
Electrical Info	220-230 V AC; 16 A; 50-60 Hz	
Dimensions & Weight	480 (W) x 900 (L) x 310 (H) mm - Approx 30kg	
Manufacturer	Lynton, UK	

Lynton Lasers Ltd

Lynton House, Manor Lane,
Holmes Chapel, Cheshire, CW4 8AF

Tel: +44 (0)1477 536 977

Fax: +44 (0)1477 536 978

info@lynton.co.uk

